





ISO9001 & ISO13485 & ISO11137

Product Name:	Petri Dish with Lock			
Product Application:	Suitable for bacterial culture. Each plate has three sets of parallel sliding clasps that			
	lock the lid and plate together without separating.			
Product Materials:	Polystyrene (PS)			
product Specifications:	Temperature range: stored at room temperature			
	Shelf life: 3 Years after date of production (ensure package is in good)			
Sterilization:	Yes.			
	Sterilized by E-beam, Sterility Assurance Level: SAL=10 ⁻⁶ . The product has been			
	irradiated and dosimetrically released based on ANSI/AAMI/ISO 11137			
Pyrogens:	Non-Pyrogenic			
BSE/TSE:	These products are deemed animal free.			
Performance Testing:	Each manufacturing lot is sampled and tested in accordance with standard operating			
	procedures.			
	Appearance inspection: pass Packaging inspection: pass			
Features	High clarity polystyrene			
	2. Flat transparent surface			
	3. Numeric mark on the bottom for easy identification			
	4. PA/PE packaging			
	5. Vent for gas exchange			
	6. Stackable for easy storage and handling			
	7. Parallel Sliding Clasps Design:			
	Each dish has three sets of parallel sliding clasp, compare with standard dish, you			
	don't need to worry about the lid or bottom falling off if you don't hold it properly.			
	8. Non-pyrogenic			
	9. Sterilized by E-Beam, SAL=10-6			

Operation approach:



Cover the lid, rotate the lid clockwise in parallel, until the clasps slide into the slots to

Head Office

No. 530, Xida Road, Meicun Industrial Park, Xinwu District, Wuxi, Jiangsu, China

Tel: +86+ 510-6800 6788 Email: info@nest-wuxi.com

Online: www.cell-nest.com



Rotate the lid counterclockwise in parallel until the clasps slip out of the slots to unlock, then the lid can be taken away.

Oversea

NEST USA (New Jersey/ Phoenix) NEST scientific 株式会社 (Yokohama,Japan) NEST Scientific Europe B.V (Netherlands) Nest Scientific (MENA) FZE (Sharjah, United Arab Emirates)








ISO9001 & ISO13485 & ISO11137

Product Range

Cat. No.	Product Description	/Pack	/Case
752101	90×15 mm Petri Dish with Lock	20	500

Technical Drawing

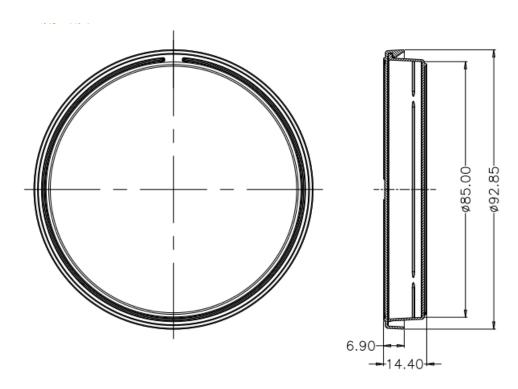

Product picture



Head Office

No. 530, Xida Road, Meicun Industrial Park, Xinwu District, Wuxi, Jiangsu, China

Tel: +86+ 510-6800 6788 Email: info@nest-wuxi.com

Online: www.cell-nest.com

Oversea

NEST USA (New Jersey/ Phoenix) NEST scientific 株式会社 (Yokohama,Japan) NEST Scientific Europe B.V (Netherlands) Nest Scientific (MENA) FZE (Sharjah, United Arab Emirates)